
# BRAUN

## Oral-B

### D 6011

### D 6511
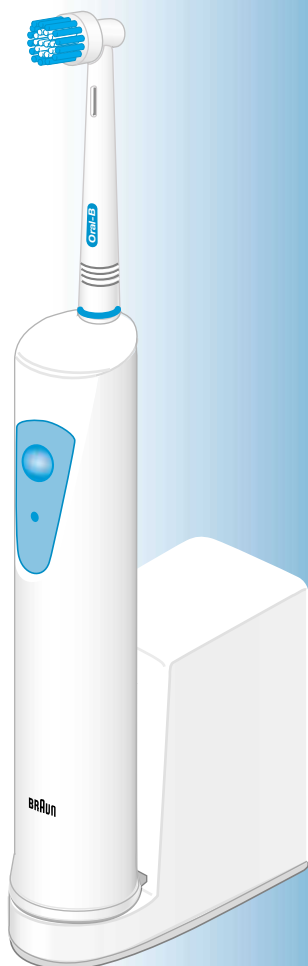
# Plak Control

### Type 4728

#### Braun Oral-B Plak Control D 6011, D 6511

Garantiekarte  
Guarantee Card  
Carte de garantie  
Carta di garanzia  
Garantiebewijs  
Garantibevis  
Köpbevis  
Takuukortti

Kaufdatum	Købsdato
Date of purchase	Kjøpsdato
Date d'achat	Inköpsdatum
Data d'acquisto	Ostopäivä
Koopdatum	

Service notes

Stempel und Unterschrift des Händlers  
Stamp and signature of dealer  
Cachet et signature du commerçant  
Timbro e firma del negozio  
Stempel en handtekening van de handelaar  
Forhandlerens stempel og underskrift  
Stempel og underskrift av forhandleren  
Återförsäljares stämpel och underskrift  
Myyjän leima ja allekirjoitus

#### Braun Oral-B Plak Control D 6011, D 6511

Registrierte Karte  
Registration Card  
Carte de contrôle  
Cartolina di registrazione  
Registratiekaart  
Registreringskort

Kaufdatum	Koopdatum
Date of purchase	Købsdato
Date d'achat	Inköpsdatum
Data d'acquisto	

Name und vollständige Anschrift des Käufers  
Name and full address of purchaser  
Nom et adresse complète de l'acheteur  
Nome e indirizzo completo dell'acquirente  
Naam en volledig adres van de koper  
Kjøperens navn og adresse  
Köparens namn och fullständiga adress

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


Deutsch	4, 13, 15, 18
English	5, 13, 15, 18
Français	6, 13, 15, 18
Italiano	7, 13, 16, 18
Nederlands	8, 13, 16, 18
Dansk	9, 13, 16, 18
Norsk	10, 14, 17, 18
Svenska	11, 14, 17, 18
Suomi	12, 14, 17, 18

## Braun Infoline

- (D)** Haben Sie Fragen zu diesem Produkt?  
Rufen Sie an: (in Deutschland und Österreich zum Nulltarif)  
**00 800 27 28 64 63**  
**00 800 BRAUNINFOLINE**
- (CH)** **08 44 -88 40 10**
- (NL)** Heeft u vragen over dit produkt?  
Bel. Braun Consumenten-infolijn:  
**(070) 4 13 16 58**
- (B)** Vous avez des questions sur ce produit?  
Appelez Braun Belgique  
**(02) 711 92 11**
- (DK)** Har du spørgsmål om produktet? Så ring  
**70 15 00 13**
- (N)** Spørsmål om dette produktet? Ring  
**88 02 55 03**
- (S)** Frågor om apparaten? Ring Kundservice  
**020 - 21 33 21**
- (FIN)** Onko Sinulla kysyttävää tuotteesta? Soita  
**0203 77877**
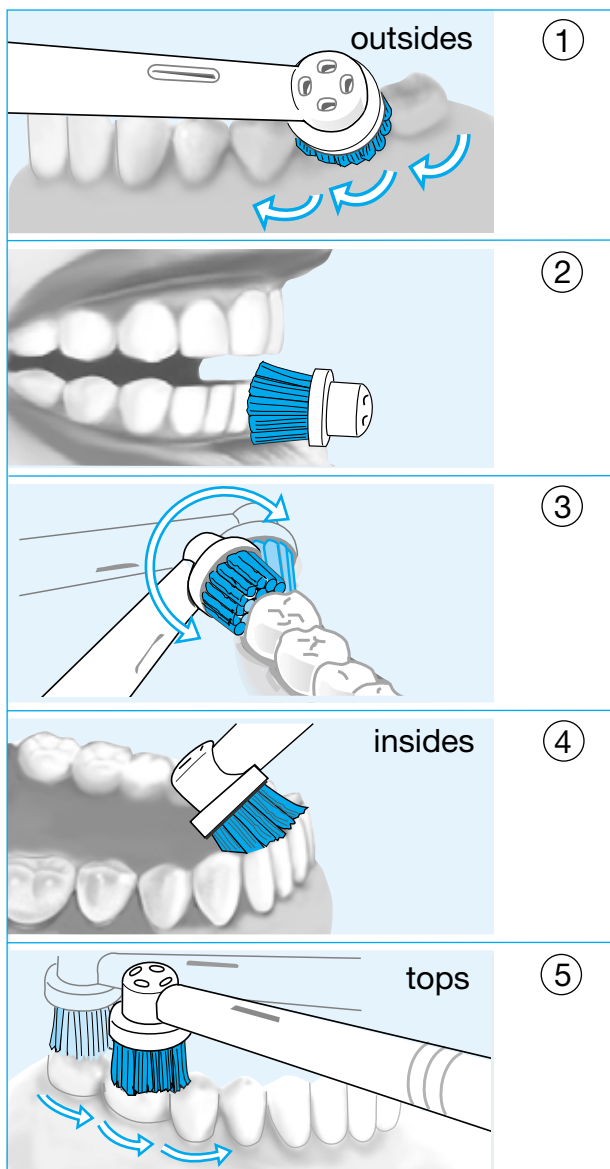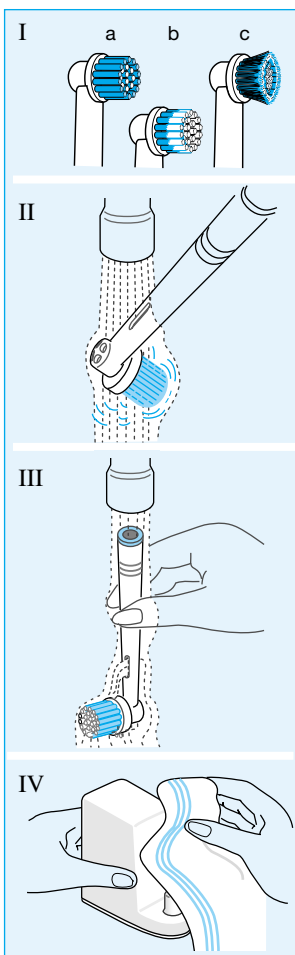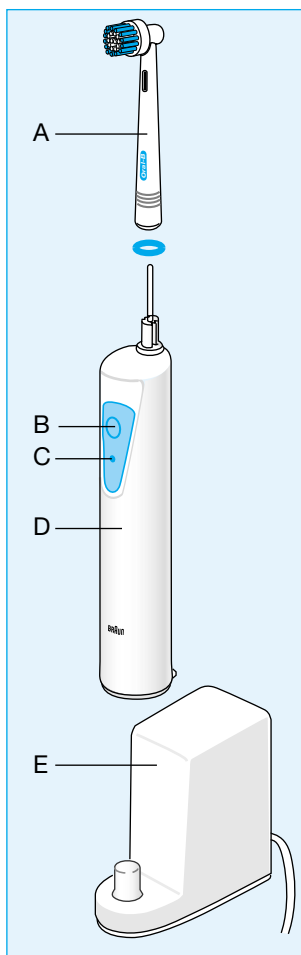
## English

The Braun Oral-B Plak Control has been developed together with leading dental professionals to effectively remove plaque from all areas of your teeth. It has been designed for daily use for the whole family.

Please read the use instructions carefully before first use.

### Important

Periodically check the cord for damage. Should this be the case, take the charging unit (E) to a Braun Service Centre. A damaged or non-functioning unit should no longer be used.

## Description

- A Brushhead
- B On switch
- C Off switch
- D Handle
- E Charging unit

## Specifications

Power supply: see bottom of charging unit  
Handle voltage: 2.4 V

## Connecting and charging

The plaque remover is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit.
- A full charge takes 16 hours and will provide an operation time of approx. 35 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the plaque remover at full power. Overcharging is impossible.

### Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

## Using the plaque remover

The plaque remover can be used with any leading brand of toothpaste. However, we do not recommend toothpastes with whitening agents.

- To avoid splashing, guide the brushhead to your teeth before switching the appliance on.
- Slowly guide the brushhead from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.
- To thoroughly remove plaque, brush for at least two minutes.

### Timer

(D 6511 only)

An incorporated timer memorises the elapsed brushing time. After two minutes the timer signals by a short stutter.

During the first days of using the plaque remover, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or hygienist.

Dentists recommend that you replace any toothbrush every 3 months, since worn bristles do not remove plaque as effectively as new ones. Therefore, the brushhead is provided with INDICATOR® bristles (I.a). With proper brushing twice a day for two minutes, the colour will disappear half-way down the bristles in a 3-month period, provided you use toothpaste (I.b).

If the bristles splay before the colour recedes (I.c), you are exerting too much pressure on teeth and gums.

### Replacement brushheads

This product can be used with all Braun Oral-B Plak Control / Plaque Remover brushheads.

## Cleaning

After use, rinse the brushhead for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brushhead. Rinse both parts separately under running water (III) and wipe them dry.

From time to time, also clean the charging unit with a damp cloth.

## Environmental notice

This product contains nickel-cadmium rechargeable batteries. In the interest of protecting the environment, do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Centre.








However, if you prefer to dispose of the battery yourself, please refer to the section «Battery removal» (page 13).

Subject to change without notice.

This product conforms to the EMC-Requirements as laid down by the Council Directive 89/336/EEC and to the Low Voltage Regulation (73/23 EEC).

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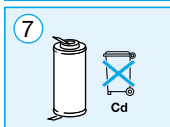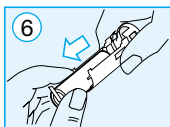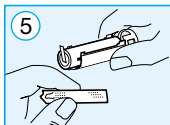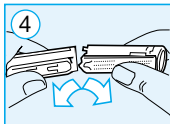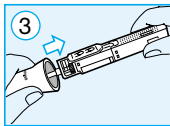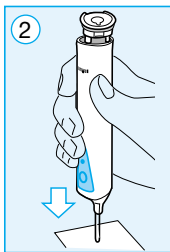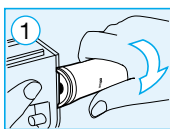## English

### Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution: Opening the handle will destroy the appliance and invalidate the guarantee.

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## English

### Guarantee

We grant a 1 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in material or workmanship either by repairing, or exchanging parts or exchanging the whole appliance as we may choose.

Not covered by this guarantee: damage arising from improper use (operation with incorrect current/voltage, plugging it into an unsuitable power point, breakage), normal wear and defects which have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

The guarantee only comes into force if the date of purchase is confirmed by the dealer's stamp and signature on the guarantee card and registration card.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor and where no import restrictions or other legal regulations hinder or prevent service being given under guarantee.

Service under the guarantee does not affect the expiry date of the guarantee. The guarantee on parts or entire appliances which are exchanged ends when the guarantee on this appliance expires.

In the event of a claim under this guarantee, hand in the complete appliance (charging unit and handle) with the valid guarantee card to your dealer, or send both to the nearest Braun Customer Service Centre.

All other claims including for damages resulting from this guarantee are excluded unless our liability is legally mandatory.

Those claims arising out of the sales contract with the vendor are not affected by this guarantee.

### For UK only:

This guarantee in no way affects your rights under statutory law.

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