### English

Our products are engineered to meet the highest standards of quality, functionality and design. We hope you thoroughly enjoy your new Braun Oral-B appliance.

The Braun Oral-B Plak Remover was developed together with leading dental professionals to effectively remove plaque from the critical areas of your teeth. Plaque is the bacterial film which forms on your teeth and can lead to tartar, gum disease and caries if not removed thoroughly.

### Caution

Check the cord occasionally for any damage. Should this be the case, unplug the charging unit ① immediately and take it to a Braun Service Centre. Your Braun Service Centre will replace the complete charging unit (free of charge within the guarantee period); the charging unit cannot be repaired for safety reasons.

A damaged or non-functioning unit may no longer be used.

### Description

- ① Charging unit
- ② Two screw openings
- ③ Colour rings
- ④ Brushheads
- ⑤ MDrive-shaft
- ⑥ Charge indicator
- ⑦ Switch
- ⑧ Rechargeable power unit
- ⑨ Roll-stop
- ⑩ Charging base

### Specifications

Charging unit: 200 - 240 volts ac (~), 50/60 Hz (without converting), 1 watt. Rechargeable power unit: 2.4 volts, maximum operation time approx. 45 minutes when fully charged.

Your Braun Oral-B Plaque Remover has inductive charging, which eliminates bothersome electrical hazards when in use. Even children will have no problem handling it. It is electrically absolutely safe, waterproof and can therefore be stored in the bathroom without hesitation.

Oscillating brush movements.  
Angle of oscillation approx. 70°.

Before using the first time, charge the power unit: if your voltage is 220 - 240 volts for at least 16 hours.

Notice: When the rechargeable cells are running empty the plaque remover automatically turns off. This protects the rechargeable cells from deep discharging. This signals that the power unit of the plaque remover should be charged.

This product contains nickel-cadmium cells. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Centre or at appropriate collection sites provided in your country.

Please, also see the insert «Cell disposal at the end of the useful life of the appliance».

### Wall fastening

Fasten wall plugs and screws in the wall. However, turn screws into wall plugs only so far that screw heads can be stuck into the two openings ② in the back of the charging unit ①. By pulling the charging unit ① down, the screws will be fastened in the two slits over the screw openings ②.

You have to pull out only as much cord from the cord storage as you need to the mains plug.

### Connecting

The charging unit ① can be connected to any voltage outlet – even to shaver sockets – between 220 - 240 volts ac (~) without converting. Place the power unit ⑧, facing the front, on the charging base ⑩.

The green charge indicator ⑥ will light up at 220 - 240 volts, whenever the power unit is inserted and recharging on the plugged-in charging unit. The power unit may warm up somewhat. This has no significance on the unit's operation. Over-charging is impossible. To maintain your plaque remover's full power capacity, the power unit should always be stored in the plugged-in charging unit after use.

**Important:** The power unit ⑧ switches off automatically when placed in the charging unit ① to permit recharging. Make sure that there is electrical power in the

socket when you recharge your power unit. Sockets in mirror cabinets over wash basins or shaver sockets in bathrooms may be without electrical power when the light is turned off.

### Switch

The power unit ⑧ is switched on and off by pressing the circular marking on the blue switch ⑦ (b). When placed in the charging unit, the power unit is automatically switched off. The blue switch is waterproof and can be cleaned under running water.

### Preserving the rechargeable batteries

To maximize the capacity and life of the rechargeable batteries, they should be fully charged (for 16 hours on 220 - 240 volts) before the first use and then fully discharged by using. Then recharge to full capacity again. This discharging/charging procedure should be done every six months.

### Use

After the power unit has been fully charged, attach one of the brushhead ④ to the power unit (a). Before doing so, you may fix one of the coloured rings ③ onto the end of the brushhead. The brushhead must face the same direction as the switch ⑦. The roll-stop ⑨ prevents the power unit ⑧ from rolling off a flat surface.

Before brushing your teeth, rinse your mouth thoroughly and wet the brushhead. You may use regular toothpaste. To avoid splashing,

switch the power unit on after guiding the brushhead to your teeth. Then, move the brushhead along the chewing surfaces of your teeth.

After that, hold the brushhead at a 90° angle and guide it along the outsides and then the insides of your teeth. Be sure to get into the spaces between your teeth as well as along the gumline. The Braun Oral-B Plak Remover will automatically provide the right cleaning action. You do not need to press hard. To thoroughly remove plaque, brush for at least 2 minutes.

After cleaning your teeth, rinse your mouth thoroughly again. During the first few days of using the Braun Oral-B Plaque Remover, your gums may bleed slightly, particularly if they normally tend to do so. After a few days, the bleeding will generally stop. If your gums continue to bleed after 2 weeks, you should see your dentist.

The brushes have soft, end-rounded bristles in order to protect your gums. Clinical tests have shown that the best cleaning effect is achieved by gently moving the brush from tooth to tooth, without exerting pressure. A change in the shape of the brush is normal after use. This is influenced by your individual cleaning technique, but does not affect the efficiency of the brush.

### INDICATOR brushheads

The INDICATOR brushhead (e / I) helps you monitor your brushing performance.

Dentists recommend that you replace any toothbrush every 3 to 4 months, since worn bristles do not clean plaque away as effectively as new ones. With proper brushing, twice a day for two minutes, the green colour will wear half-way down the bristles in this 3 to 4 month period, provided you use toothpaste (e / II).

If the bristles begin to splay before the green colour starts to recede it means that you are exerting too much pressure on your teeth and gums, and wearing out the brushhead too quickly (e / III).

The INDICATOR brushhead helps to improve your brushing performance, thus optimizing exceptional plaque removal and the life of the brushhead.

Replacement brushheads are available at your dealer or Braun Service Centres.

### Cleaning

To maintain the function of the brushhead ④ rinse thoroughly for several seconds with the power unit ⑧ switched on under running water (not warmer than 80 °C) after use (c).

Clean the brushheads and the power unit by simply holding them under running water (d). Wipe the power unit dry. After cleaning place the power unit ⑧ back on the charging base ⑩. The green charge indicator ⑥ will light up at 220 - 240 volts.

Additionally, we recommend soaking the brushheads in a denture cleaner (commercially available) every 2 to 4 weeks to keep them hygienically clean. Proceed according to the manufacturer's use instructions.

Subject to change without notice.

**CE** This product conforms to the EMC-Requirements as laid down by the Council Directive 89/336/EEC.

#### **Guarantee**

We grant a 1 year guarantee on the product commencing on the date of purchase.

Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in material or workmanship either by repairing, or exchanging parts or exchanging the whole appliance as we may choose.

Not covered by this guarantee: damage arising from improper use or storage (operation with incorrect current/voltage, plugging it into an unsuitable power point, breakage), normal wear and defects which have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are under-taken by unauthorized persons and if original Braun parts are not used. The

guarantee only comes into force if the date of purchase is confirmed by the dealer's stamp and signature on the guarantee card and registration card.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor and where no import restrictions or other legal regulations hinder or prevent service being given under guarantee.

Service under the guarantee does not affect the expiry date of the guarantee. The guarantee on parts or entire appliances which are exchanged ends when the guarantee on this appliance expires.



In the event of a claim under this guarantee, hand in the complete appliance (recharger and rechargeable power unit) with the valid guarantee card to your dealer, or send both to the nearest Braun Customer Service Centre.

All other claims including for damages resulting from this guarantee are excluded unless our liability is legally mandatory.

Those claims arising out of the sales contract with the vendor are not affected by this guarantee.

#### **For UK only:**

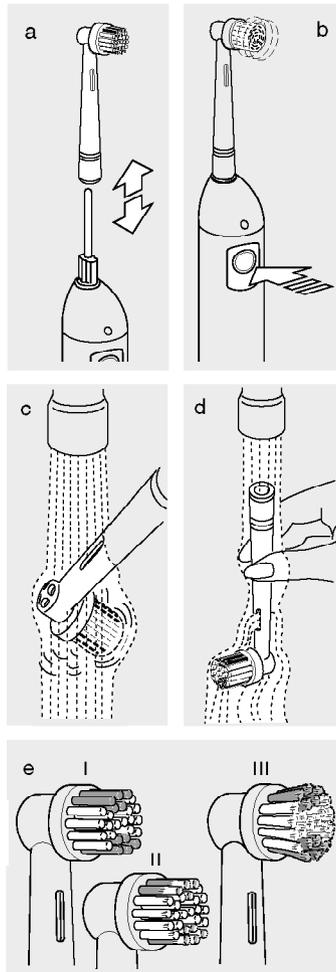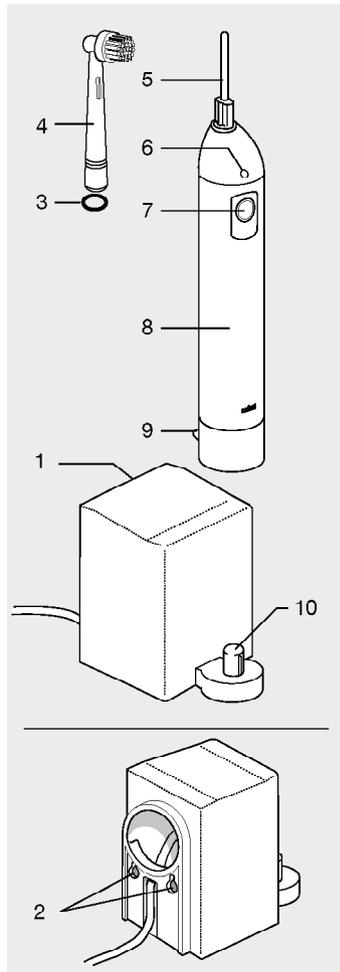This guarantee in no way affects your rights under statutory law.

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