Oral-B®

Professional Care

D 15525

D 15513

D 15513 T

D 15511

powered by **BRAUN**

Type 4729

Type 4736

Deutsch	3, 19, 22
English	4, 19, 22
Français	5, 19, 22
Italiano	6, 19, 22
Nederlands	7, 20, 22
Dansk	8, 20, 22
Norsk	10, 20, 22
Svenska	11, 20, 22
Suomi	12, 21, 22
Ελληνικά	15, 21, 22


- D** Haben Sie Fragen zu diesem Produkt?
- A** Rufen Sie an (gebührenfrei): **00 800 27 28 64 63**
- CH** **08 44 - 88 40 10**
- GB** If you should require further assistance, please call our Consumer Service line. **0 800 73 11 792** (free call)
- IRL** Helpline **1 800 509 448** (free call)
- F** Appelez le Service Consommateurs Oral-B: **0 810 365 855** (numéro azur)
- B** **02-711 92 11**
- LU** Vous avez des questions sur ce produit? Appelez: **4 00 50 51**
- I** Servizio consumatori: **02-667 86 23**
- NL** Heeft u vragen over dit produkt? Bel Braun/Oral-B Consumenten- infolijn: **0800-4 45 53 88**
- DK** Har du spørgsmål om produktet? Så ring **70 15 00 13**
- N** Spørsmål om dette produktet? Ring **22-72-88-10**
- S** Frågor om apparaten? Ring Kundservice **020 - 21 33 21**
- FIN** Onko Sinulla kysyttävää tuotteesta? Soita **09 452 872**
- GR** Τηλέφωνο για πληροφορίες **210 9478700**


Internet:
www.oralb.com
www.braun.com

3-729-162/00/III-05/M
D/GB/F/I/NL/DK/N/S/FIN/GR
Printed in Germany

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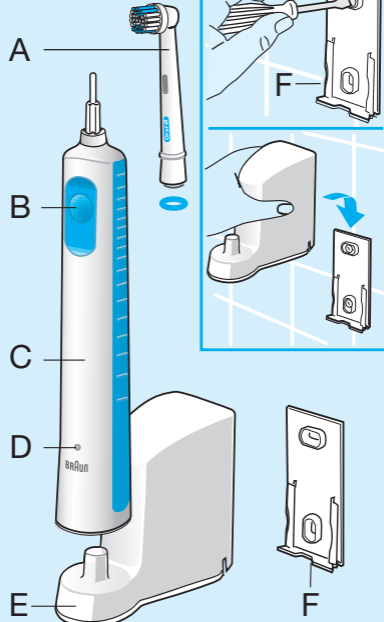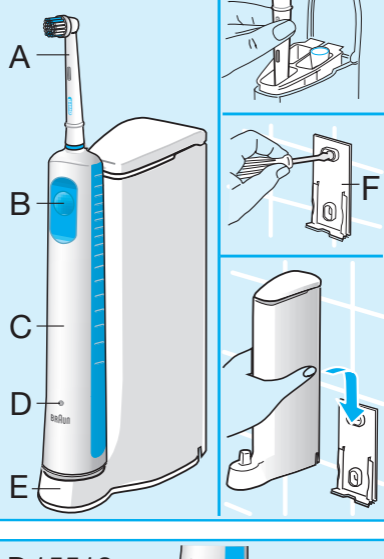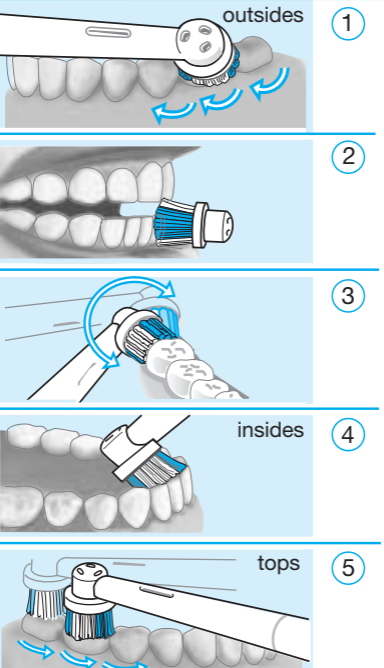
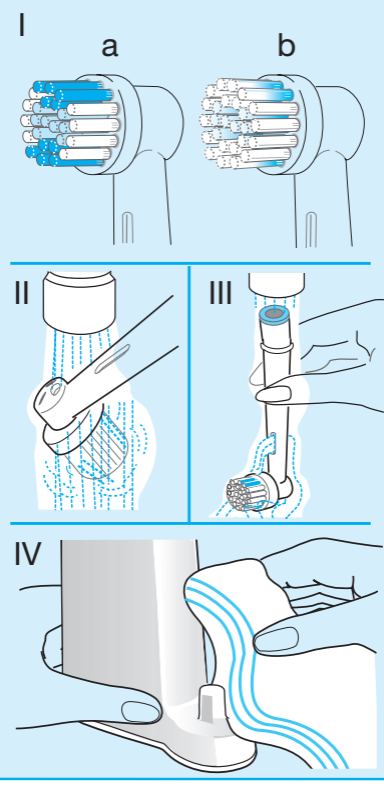
English

Please read the use instructions carefully before first use.

Important

Periodically check the cord for damage. Should this be the case, take the charging unit (E) to a Braun Service Centre. A damaged or non-functioning unit should no longer be used.

Description

- A Brushhead
- B On/off switch
- C Handle
- D Charge indicator
- E Charging unit
- F Wall holder (D 15511 and D 15525 only)

Specifications

Power supply: see bottom of charging unit

Handle voltage: 1.2 V

Connecting and charging

The toothbrush is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit.
- A full charge takes approx. 12 hours and will provide an operation time of approx. 30 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the plaque remover at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Using the toothbrush

The toothbrush can be used with any leading brand of toothpaste.

- To avoid splashing, guide the brushhead to your teeth before switching the appliance on.
- Slowly guide the brushhead from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.
- To thoroughly remove plaque, brush for at least two minutes.

Memory timer

An incorporated timer memorises the elapsed brushing time, even when the handle is briefly (max. 30 sec) switched off. After two minutes, the timer signals by a short stutter.

Pressure Sensor

If you apply too much pressure, the oscillating movement will continue, but the pulsation will stop. There is a noticeable change in feeling and sound when the Pressure Sensor is activated.

During the first days of using the appliance, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or hygienist.

Dentists recommend that you replace any toothbrush every 3 months.

The INDICATOR[®] bristles help you monitor and improve your brushing performance. With proper brushing, twice a day for two minutes, the colour on the INDICATOR[®] bristles will disappear half-way down within 3 months, signaling it is time to change the brushhead.

Replacement brushheads

Available at your retailer or Braun Service Centres.

Cleaning

After use, rinse the brushhead for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brushhead. Rinse both parts separately under running water (III) and wipe them dry.

From time to time, also clean the charging unit with a damp cloth.

Environmental notice

This product is provided with a nickel-hydrate battery which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Centre.

However, if you prefer to dispose of the battery yourself, please refer to the chapter «Battery removal» on page 17.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

At the end of the product's useful life, please dispose of it at appropriate collection points provided in your country.

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