
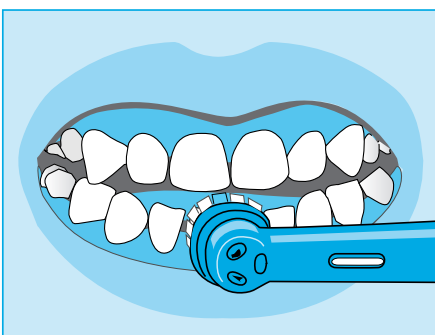
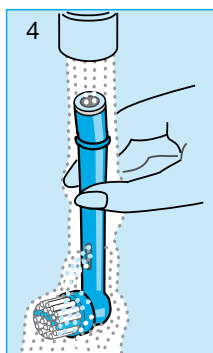
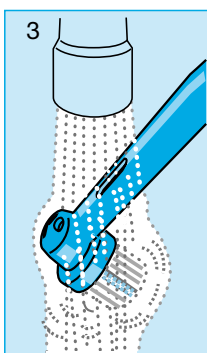
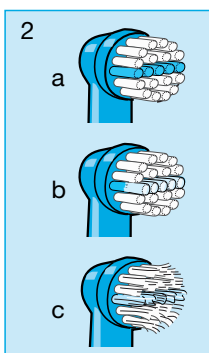
Oral-B®**AdvancePower
D 10 Kids**powered by **BRAUN****Type 4733****Braun Infoline**

Deutsch	3, 16, 19
English	4, 16, 19
Français	5, 16, 19
Español	6, 16, 19
Português	7, 16, 19
Italiano	8, 17, 19
Nederlands	9, 17, 19
Dansk	10, 17, 19
Norsk	11, 17, 19
Svenska	12, 18, 19
Suomi	14, 18, 19

Internet:www.oralb.com
www.braun.com


- (D)** Haben Sie Fragen zu diesem Produkt?
- (A)** Rufen Sie an (gebührenfrei):
00 800 27 28 64 63
00 800 BRAUNINFOLINE
- (CH)** **08 44 - 88 40 10**
- (E)** Servicio al consumidor para España:
901 11 61 84
- (P)** Serviço ao Consumidor para Portugal:
808 20 00 33
- (I)** Servizio consumatori:
(02) 6 67 86 23
- (NL)** Heeft u vragen over dit product? Bel Braun Consumenten-infolijn:
(070) 4 13 16 58
- (B)** Vous avez des questions sur ce produit? Appelez Braun Belgique
(02) 711 92 11
- (DK)** Har du spørgsmål om produktet? Så ring
70 15 00 13
- (N)** Spørsmål om dette produktet? Ring
88 02 55 03
- (S)** Frågor om apparaten? Ring Kundenservice
020 - 21 33 21
- (FIN)** Onko Sinulla kysyttävää tuotteesta? Soita
0203 77877

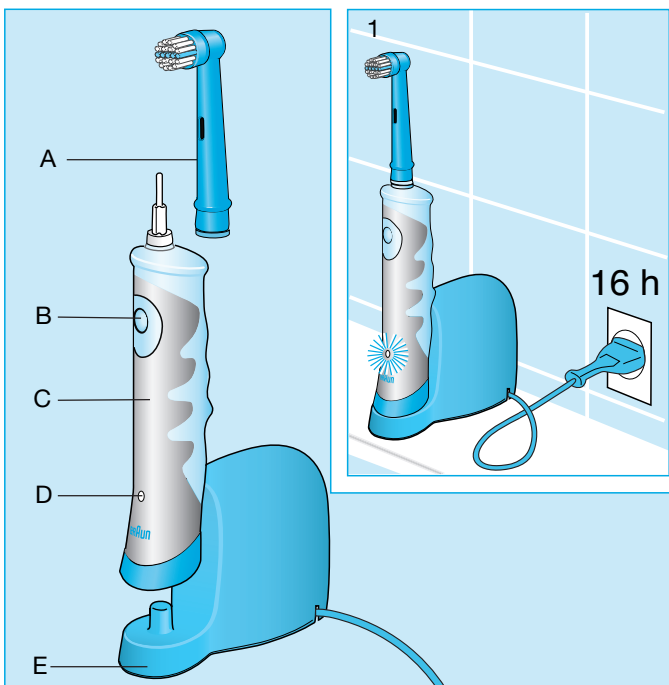outsides

insides

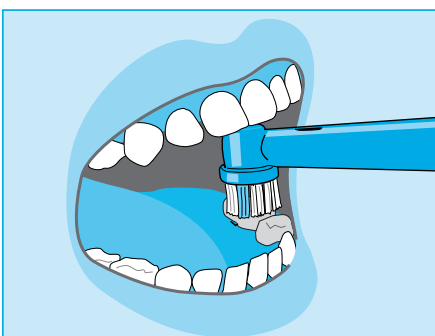tops

English

Dear parents,
The Oral-B AdvancePower Kids has been developed with leading dental professionals to effectively remove plaque from children's teeth. It will make it easier for your child to brush properly and thoroughly. We hope your child will enjoy his/her new Oral-B appliance.

This appliance is suitable for children 3 years or older. Please read the use instructions carefully and completely before showing your child how to use the appliance. We recommend that parents supervise their children when they start using the appliance in order to help them become familiar with it.

Important

Periodically check the cord for damage. Should this be the case, take the charging unit (E) to a Braun Service Centre. A damaged or non-functioning unit should no longer be used.

Description

- A Brushhead
- B On/off switch
- C Squish Grip handle
- D Charge indicator
- E Charging unit

Specifications

Power supply: see bottom of the charging unit

Handle voltage: 1.2 V

Connecting and charging

The plaque remover is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit. The charge indicator will light up (1).
- A full charge takes 16 hours and will provide an operation time of approximately one hour.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the appliance at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Removing plaque with a music timer

Use the following steps to show and explain to your child how to use the appliance properly:

- When the handle is taken from the plugged-in charging unit, a short melody will sound.
- Put a small amount of toothpaste onto the brushhead. (We recommend fluoride toothpaste.)
- First guide the brushhead to the teeth, only then turn on the handle with switch (B).

- Start with the lower teeth: Slowly guide the brushhead from tooth to tooth;
 - first the outsides,
 - then the insides,
 - finally the chewing surfaces.

Do not scrub or press too hard. The Oral-B AdvancePower Kids automatically provides the proper cleaning action.

- After one minute of operation, the brush movement will be interrupted by a short melody. This tells your child to move on to the upper teeth. Clean them – again tooth by tooth – until the next melody plays (after another minute of operation). This way you can make sure that your child brushes his/her teeth for at least 2 minutes as recommended by dental professionals.
- Turn off the handle with switch (B) while the brushhead is still in the mouth.

We recommend that you replace the brushhead every three months. The brushhead is provided with a row of blue INDICATOR® bristles (2a). With proper brushing twice a day for two minutes, the colour will disappear half-way down these bristles in a 3-month period (2b), provided toothpaste is used. If bristles begin to splay before the colour starts to recede, this indicates that your child exerts too much pressure on teeth and gums (2c).

Replacement brushheads EB10

Available at your retailer or at Braun Service Centres.

Cleaning

After use, rinse the brushhead for several seconds under running water with the handle switched on (3). Then switch off the handle and take off the brushhead. Rinse both parts separately under running water (4) and wipe them dry before putting the handle back onto the charging unit.

Environmental notice

This product is provided with a nickel-hydrate rechargeable battery which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Centre.


However, if you prefer to dispose of the battery yourself, please refer to the chapter «Battery removal» on page 15.

Subject to change without notice.

This product conforms to the EMC-Directive 89/336/EEC and to the Low Voltage Regulation (73/23 EEC).

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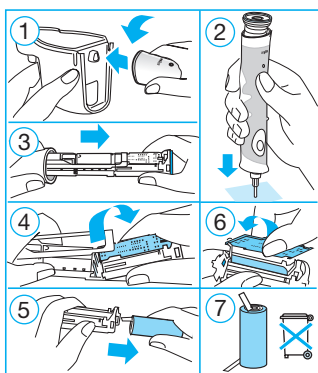English

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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